



ULTRA-THIN SPECIALIST CONTACT LENS FOR KERATOCONUS AND IRREGULAR CORNEAS



INCREASED OXYGEN PERMEABILITY

Helps maintain a healthy cornea – ideal for post-graft cases and reverse geometry corneas.



MAXIMUM COMFORT

0.20 CT improves comfort allowing the lens to stay on-eye for longer.



IMPROVED CORNEAL DRAPAGE

Enhances the drape effect in the periphery which improves fit for complex cases such as nipple cones and PMD.

PARAMETERS

Base Curves	7.40 mm to 9.40 mm (0.20 steps)
Diameters	14.50 mm (0.50mm steps)*
Spherical Powers	+/- 30.00DS (0.25 steps)
Cylinder	-0.50D to -15.00D
Axes	1° to 180°
Material	

¹ Efofilcon A, 74% Water**3 month replacement

² Filcon II 3, 77% Water 12 month replacement

* Diameters of 14.00mm, 15.00mm and 15.50mm can be manufactured to order.


** Definitive™

THE NEWEST INNOVATION FOR KERATOCONUS AND IRREGULAR CORNEAS

KeraSoft® Thin, the newest innovative addition to UltraVision CLPL's highly successful specialist range of KeraSoft® contact lenses, is designed to fit all irregular corneas, including post-refractive surgery, Pellucid Marginal Degeneration (PMD) and keratoconic conditions. The main ocular health benefit of KeraSoft® Thin is a reduced centre thickness (0.20CT) which increases oxygen permeability and improves overall comfort.

KERASOFTLENS.COM/TRAINING
+44 (0) 1525 381112



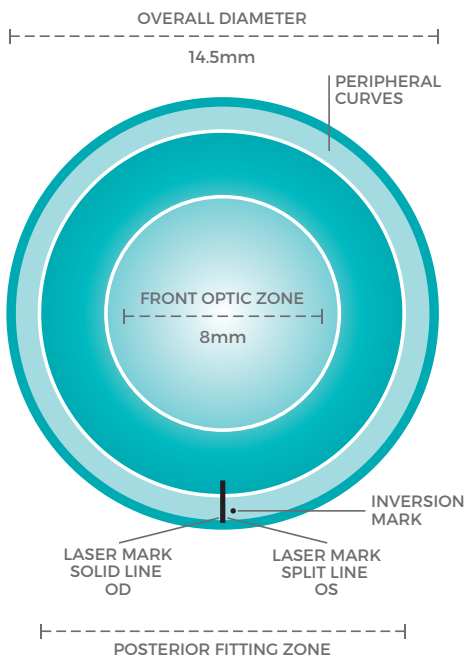

GUARANTEED FIT PROGRAM

KeraSoft® Thin contact lenses are available in a single vial with a warranty that provides for two exchanges. After successfully fitting a patient using a single lens (with warranty), you can save money by ordering replacement non-warranted lenses in a 2-pack. KeraSoft® Thin is available in the same parameter range as KeraSoft® IC and can be ordered from trials performed with the 0.40 Fitting Set.

MC036.06.16.REV.3

STANDARD FITTING SET PARAMETERS*

BASE CURVE	DIAMETER	PERIPHERY	POWER
7.80mm	14.5mm	STD	Plano
8.00mm	14.5mm	STD	Plano
8.20mm	14.5mm	STD	Plano
8.40mm	14.5mm	STD	Plano
8.60mm	14.5mm	STD	Plano
8.80mm	14.5mm	STD	Plano
8.20mm	14.5mm	FLT2	Plano
8.60mm	14.5mm	STP2	Plano

*Fitting set supplied in 77% material

ORDERING AND RETURN PROCESS

To place your order, please call any of the numbers listed below.

Warranted KeraSoft® Thin lenses being returned for credit or exchange must meet the following guidelines:

- Returned in original vial with original label
- Returned within the warranty period
- Returned with a copy of the returns form and reason for return

UltraVision CLPL

Commerce Way, Leighton Buzzard,
LU7 4RW, United Kingdom
Tel: +44 (0) 1525 381112
Fax: +44 (0) 1525 370091
Email: info@ultravision.co.uk